# nevir®

MANUAL DE INSTRUCTION

**BÁSCULA DE ANÁLISIS CORPORAL**

**NVR-3330 BLB**

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## DEAR CUSTOMER:

Welcome to use our Bluetooth Body Fat Scale.

This product is specially designed to test body fat/hydration/muscle/bone and help keep a close eye on the daily change of your body fat, hydration ,muscle, bone and weight, providing information for your reference. Based on the analysis of your height,weight and age etc., it tells your daily energy requirement and helps control your weight. Please read this instruction carefully before use.

## ABOUT SENSSUN

Bluetooth scale helps you to store and keep track of your weight on mobile devices. This scale can transmit your weight data to your smart devices via Bluetooth technology as well as make graph of your weight results, from which you are able to keep track of the trend anytime. Download SENSSUN Scales software to use Bluetooth scale and you will start your first step of weight control.

## USER-FRIENDLY TIPS

### 1. During Measuring

- A. Please always weigh/ use the scale on a solid and flat surface. The device is designed to allow auto step on. Please always ignore the first reading and only start taking effective record from the second time weighing.
- B. Always remove your shoes and socks and wipe clean your feet before using.
- C. Step on the weighing platform gently.
- D. It is advisable to take measurement at the same time of the day.
- E. Result could be misleading after intensive exercise, excessive dieting or under extreme dehydration condition.
- F. Keep away from strong electro-magnetic field when using the scale.
- G. Children under 10 or adult over 100 can only use the device under normal weighing mode.
- H. The measured data of the following people may have deviation:
  - Adult over 70
  - Body building athletes or other occupational athletes.
- I. People with the following conditions should not use the device: Pregnant women, Those with symptom of edema, or on dialysis treatment.
- J. The Fat% Hydration% content evaluated by the Fat Monitor Scale is only for reference (not for medical purpose). If your Fat% Hydration% content is over/below normal level, please consult your doctor for more advice.
- K. If failed to turn on the scale, check whether battery is installed or battery power is low. Insert or replace new battery.
- L. Do not use it for commercial purpose.

## 2. Safety Warning

- A. It could be slippery if you weigh on a wet surface. Moreover, if the weighing platform surface is glazed, never step on the edge of the scale and do keep balance during your weighing.
- B. Those who use heart pacemaker or with other implanted medical device should not use body fat testing function.

## FEATURES

- Using high precision strain gauge sensor
- Auto on function.
- Low battery and overload indication.
- Directly check personal weight, body fat, body hydration, body muscle ,d body bone percentage and suggest calorie intake requirement.
- Set password to protect personal information and store multi-suer data.
- Bluetooth wireless technology.
- Connect with different smart devices respectively.

## BATTERY INSTALLATION


Name	Illustration	Installation
2x1.5V AAA or 4x1.5V AAA Battery		AAA battery installation Please install battery according to polarity.
1x3V or 2x3V lithium battery(CR2032)		A. Take out the used battery with the aid of a cuspidal object B. Install new battery by putting one side of the battery beneath the battery baffle and then pressing down the other side. Please take note to the polarity.

Note: Please use corresponding battery. Please remove battery if it is not in use in a long time(more than three months).

## WEIGHT AND BODY FAT MEASUREMENT WITH APP


Make sure your feet on electrode or it can not get the data correctly,  
Please take off socks to ensure accuracy of data.

1. Open App. Step on scale to turn it on. Make sure the Bluetooth is on.
2. Step on scale. Weight is showed on both scale and App. After several seconds, you can get your body fat, hydration, muscle, bone and kcal on App.
3. Data will be saved automatically in history record.

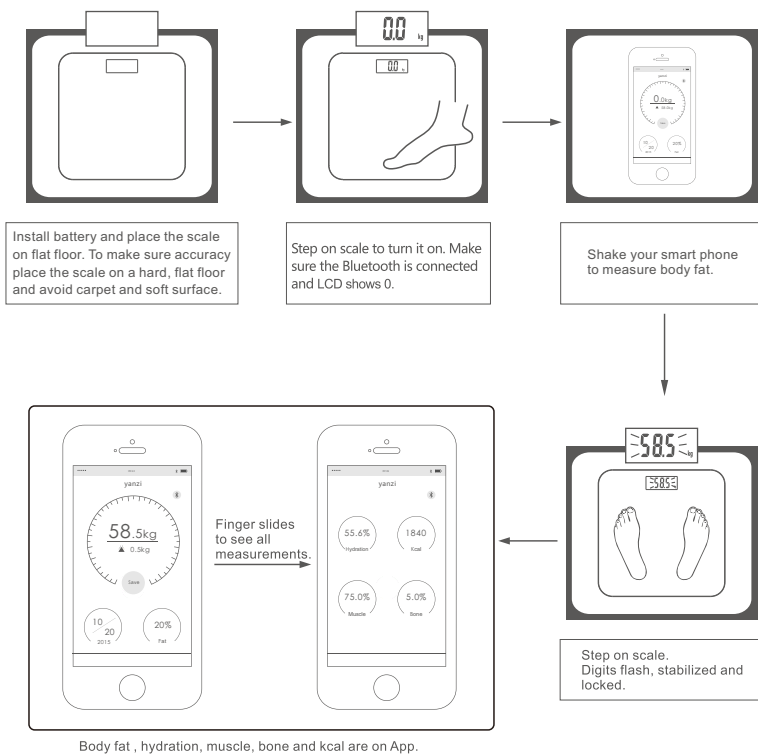
## NORMAL WEIGHING MODE

After installing battery, place the scale on flat floor, Keep stand still on scale when teting. Then wait for seconds, display shows the same weingt reading for approximate a few seconds.

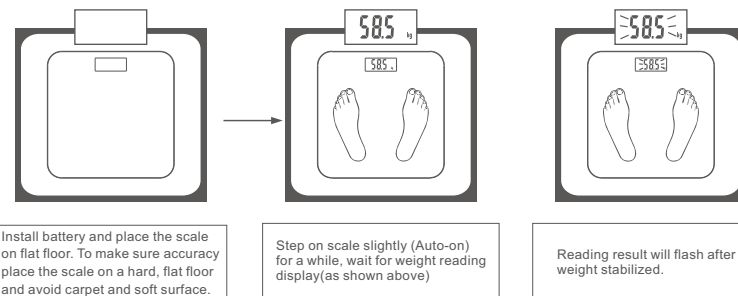
## TO CONVERT WEIGHT UNIT

If the scale is with both measuring system of Metric(kilograms) and Imperial (pounds, stones), you can convert weight unit as follow:

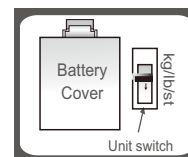
- Condition 1: if the unit switch is on the back of the scale,
- A. Find unit switch on the back of scale.
  - B. Pull up and down the unit switch to convert scale unit between kg-lb-st.

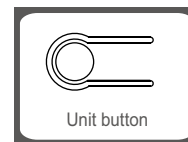- Condition 2: if the unit switch is on the back of the scale,
- A. Find unit switch on the back of scale.
  - B. Press the unit conversion button when scale is power on, LCD shows current weight unit. Press the button again to convert weight unit.

## AUTOMATICALLY TURN OFF

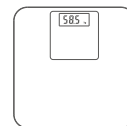Data stabilized and locked for several seconds. It will turn off automatically with no using continuously.

## ■ INDICATION

Warning indication	Name	Instructions and methods
Err	OVERLOAD INDICATION	Scale overload the max weight (Max weight is indicated on the scale.) Please step off to avoid damage.
Lo	LOW BATTERY INDICATION	Battery power is low. Please replace new battery.
C	REMEASURE	There is error. Please remeasure.

## ■ ADVISE FOR USE & CARE

1. DO NOT use the scale in humid, overheat or overcold environment. (Temperature range from 5°C to 35°C)
2. DO NOT use the scale on soft mat, lawn, slippery floor or uneven floor.
3. DO NOT drop, shock or strike the scale as it is a precise instrument.
4. Keep the scale clean. DO NOT use acid cleaning agent.

## ■ APP INSTRUCTION

1. Get the application
  - a. Scan the QR code on the right side to obtain the App.
  - b. At the "App Store" or "Google Play", search keyword "SenssunHealth" to download our application.
2. Bluetooth settings
  - a. Make sure that the Bluetooth is turned on before use
  - b. If not enabled, please go to the "Setting < Bluetooth" menu to turn it on.
  - c. Press the scale platform to turn on.
3. Create a new user
  - a. Click "SenssunHealth" icon, you will enter the "user settings".
  - b. Complete the information to enter weighing interface.
4. How To Use The Scale With Application
  - a. Press scale to turn it on
  - b. Open App and make sure the "  " on main interface is blue, which means scales connect the App.
  - c. Shake the smart phone to transfer user's information to scale.
  - d. Step on scale.
  - e. Weight will be showed on the scale and App. After several seconds, body fat, hydration, muscles, bone percentage, and kcal suggestion will be showed following upon weight measurement.
  - f. Click "Save" to keep data in history.
  - g. You can review the history, add new users and set reminders for your weighing habit.

### Warm Tips:

Requirements for mobile devices:

Bluetooth version: Bluetooth 4.0 or above

Apple System: iPhone 4S or later, iOS 6.0 or later

Android System: Android 4.3 or later

## ■ FAQ

The software can't obtain weighing results from Bluetooth scale.

- Please confirm the Bluetooth icon on the software interface is in dark blue, showing connected successfully.
- Confirm the Bluetooth is on in your device.
- Confirm the Bluetooth scale is on within the valid range of signal coverage.
- Relaunch the application.
- Turn off the Bluetooth and on again, then bridge the Bluetooth connection between scale and smart phone.

Fail to connect to your device

- Please make sure only one device is setting up connection with the Bluetooth scale.
- Please make sure your device is equipped with Bluetooth 4.0 or above.

LCD has no responses

- Please first check the batteries and replace them if necessary.
- Repairs may only be carried out by our authorized customer service or dealers.

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Importador: **NEVIR, S.A.**
Dirección del importador: NEVIR, S.A.